


QUICK START GUIDE

BEFORE YOU BEGIN!

- Scan the QR code with your mobile phone or visit accustarlabs.com/water to register device and fill out Data Form.
- Ensure that testing guidelines are followed. Scan the QR code for more information.

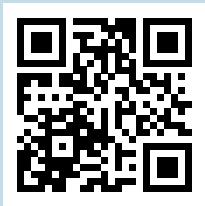

• **RECORD YOUR DEVICE NUMBER(S):** _____

KEEP FOR YOUR RECORD

Start Test

PRIOR TO EXPIRATION DATE

- 1 Test Instructions**
Read the instructions thoroughly and scan the QR code or visit our website for more detailed instructions and supplementary video.



- 2 Remove any aeration devices or filters from the tap.**
Run the cold water until fresh water is flowing from the well.

- 3 Slowly fill a bowl with the spigot underwater.**



- 4 Submerge the vial and cap open side up until they fill with water.**



While still underwater, screw the cap on the vial.



QUESTIONS? CALL US!

800-767-3703

You must have your device number(s) and zip code of the property tested.

Details and videos online:

 accustarlabs.com/water

End Test

- 5 Remove the vial from the water and check for bubbles.**
If there are any bubbles, empty the vial and repeat the process. Note: If conducting a duplicate test, immediately take the second sample from the same spigot.



Place vial back into the return box. Double kits will have two vials.

- 6 Return Device to Laboratory**
Using FedEx or UPS (USPS not recommended), return the device to the lab. Ensure that your device reaches the lab within **4 days** of the sample date.



Getting Your Report

You will receive your report via email the next business day after your device arrives at the lab.

You may also view your report at accustarlabs.com/water, where you will find additional information and resources regarding your report.

TIME SENSITIVE!

UPON COMPLETION, RETURN
TEST AS SOON AS POSSIBLE.
SEE REVERSE SIDE FOR DETAILS.

